

Rat Matrix metalloproteinase-2 ELISA Kit (2 plates)

Catalogue No.: BEK-2187-2P

Description: The rat Matrix metalloproteinase-2 (MMP-2) Kit is a sandwich ELISA. The capture antibody is a monoclonal rat MMP-2 antibody pre-coated onto the 96-well strip plates provided in the kit. Test samples and standards of known MMP-2 concentration are added to these wells and allowed to complex with the bound MMP-2 antibody. A biotinylated anti-rat MMP-2 polyclonal antibody is then added. This detection antibody binds to the antigen thus completing the sandwich. After washing, an enzyme Avidin-Biotin-Peroxidase complex (ABC) is added which binds to the second antibody. The peroxidase substrate TMB is added to induce a coloured reaction product. The intensity of this coloured product is directly proportional to the concentration of MMP-2 present in the samples. The purpose of this kit is the in-vitro quantitative determination of rat MMP-2 in samples such as sera, plasma, and cell culture supernates. This kit has been configured for research use only and is not to be used in diagnostic or clinical procedures.

Batch No.: See product labels

Antigen: The matrix metalloproteinases (MMPs) are a large family of zinc endopeptidases. All MMPs are synthesized as inactive proenzymes. The activation of these proenzymes is a critical step that leads to degradation of extracellular matrix components such as fibronectin and collagen type III. MMP2 is activated by other MMPs such as MMP14 and MMP16. MMP2 is known to degrade collagen types IV, V, VII and X as well as gelatin type I.

Other Names: 72 kDa type IV collagenase; 72 kDa gelatinase; Matrix metalloproteinase-2; MMP-2; Gelatinase A; TBE-1; MMP2; CLG4A;

Accession: P33436 MMP2_RAT;

Specificity: MMP-2

Storage: Store at 4°C

Kit components: The ELISA kit box contains 2 x 96 pre-coated strip plates, protein standards, detection reagents, substrate buffer and precise instructions.

Range: 625 pg/ml - 40,000 pg/ml


Sensitivity: < 10pg/ml

Kit protocol: Please refer to our online product listing for current protocol/MSDS versions.

IX: Typical Standard Curve (for reference only, not to be used for actual data)

| Concentration pg/ml | 0.0 pg/ml | 625 pg/ml | 1250 pg/ml | 2500 pg/ml | 5000 pg/ml | 10000 pg/ml | 20000 pg/ml | 40000 pg/ml |
|------------------------|-----------|-----------|------------|------------|------------|-------------|-------------|-------------|
| O.D. | 0.086 | 0.192 | 0.233 | 0.335 | 0.524 | 0.833 | 1.478 | 2.119 |

This standard curve is for demonstration purposes only. A standard curve should be generated for each assay.

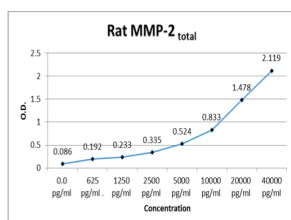

FOR RESEARCH USE ONLY