



Mouse Osteopontin ELISA Kit (2 plates)

Catalogue No.: BEK-2084-2P

Description: The mouse Osteopontin (OPN) Kit is a sandwich ELISA. The capture antibody is a polyclonal mouse OPN antibody pre-coated onto the 96-well strip plates provided in the kit. Mouse test samples and standards of known OPN concentration are added to these wells and allowed to complex with the bound OPN antibody. A biotinylated mouse OPN polyclonal antibody is then added. This detection antibody binds to the antigen thus completing the sandwich. After washing, an enzyme Avidin-Biotin-Peroxidase complex (ABC) is added which binds to the second antibody. The peroxidase substrate TMB is added to induce a coloured reaction product. The intensity of this coloured product is directly proportional to the concentration of OPN present in the samples. The purpose of this kit is the in-vitro quantitative determination of mouse OPN in samples such as sera, plasma, tissue lysates and cell culture supernates. This kit has been configured for research use only and is not to be used in diagnostic or clinical procedures.

Batch No.: See product labels

Other Names: OPN; 2AR; Bone sialoprotein 1; Calcium oxalate crystal growth inhibitor protein; Early T-lymphocyte activation 1 protein; Minopontin; Secreted phosphoprotein 1; SPP-1; Eta-1; Op; Spp1;

Accession: P10923 OSTP_MOUSE;

Specificity: Mouse OPN

Storage: Store at 4°C

Kit components: The ELISA kit box contains 2 x 96 pre-coated strip plates, protein standards, detection reagents, substrate buffer and precise instructions.


Range: 156 pg/ml - 10,000 pg/ml

Sensitivity: < 10 pg/ml

Kit protocol: Please refer to our online product listing for current protocol/MSDS versions.

IX: Typical Standard Curve (for reference only, not to be used for actual data)

Concentration pg/ml	0.0	156	312	625	1250	2500	5000	10000
O.D.	0.002	0.082	0.160	0.304	0.560	1.030	1.453	1.898

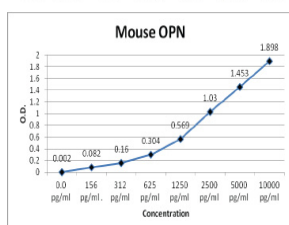

This standard curve is for demonstration purposes only. A standard curve should be generated for each assay.

FOR RESEARCH USE ONLY