

## Human Interleukin 17 (IL-17) ELISA Kit (2 plates)

**Catalogue No.:** BEK-2053-2P

**Description:** The human Interleukin 17 (IL-17) Kit is a sandwich ELISA. The capture antibody is a polyclonal human IL-17 antibody pre-coated onto the 96-well strip plates provided in the kit. Human test samples and standards of known IL-17 concentration are added to these wells and allowed to complex with the bound IL-17 antibody. A biotinylated human IL-17 polyclonal antibody is then added. This detection antibody binds to the antigen thus completing the sandwich. After washing, an enzyme Avidin-Biotin-Peroxidase complex (ABC) is added which binds to the second antibody. The peroxidase substrate TMB is added to induce a coloured reaction product. The intensity of this coloured product is directly proportional to the concentration of IL-17 present in the samples. The purpose of this kit is the in-vitro quantitative determination of human IL-17 in samples such as sera, plasma, tissue lysates and cell culture supernates. This kit has been configured for research use only and is not to be used in diagnostic or clinical procedures.

**Batch No.:** See product labels

**Other Names:** Interleukin-17; IL-17; IL17; Interleukin-17A; Cytotoxic T-lymphocyte-associated antigen 8; CTLA-8; IL17A; IL17;

**Accession:** Q16552 IL17\_HUMAN;

**Specificity:** Human IL-17

**Storage:** Store at 2-8C

**Kit components:** The ELISA kit box contains 2 x 96-well pre-coated strip plates, protein standards, detection reagents, substrate buffer and detailed protocols.


**Range:** 31.5 pg/mL - 2000 pg/mL

**Sensitivity:** < 1 pg/ml

**Kit protocol:** Please refer to our online product listing for current protocol/MSDS versions.

IX: Typical Standard Curve (for reference only, not to be used for actual data)

Concentration	0.0pg/ml	31.2pg/ml	62.5pg/ml	125pg/ml	250pg/ml	500pg/ml	1000pg/ml	2000pg/ml
O.D.	0.030	0.082	0.137	0.219	0.440	0.847	1.407	2.242

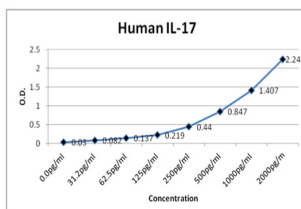

This standard curve is for demonstration purposes only. A standard curve should be generated for each assay.

FOR RESEARCH USE ONLY